KLORİT

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REFERANSLAR

1. Tietz, N.W., Fundamentals of Clinical Chemistry, p. 940, W.B. Saunders Co., Philadelphia, 1987.
2. Tietz NW. Clinical Guide to Laboratory Test. 2nd ed. Philadelphia, PA: WB Saunders Company; 1995,52.
3. Tietz NW. Clinical Guide to Laboratory Tests. 3rd ed. Philadelphia, PA: WB Saunders Company; 1995:88-91.

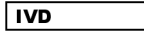
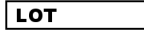
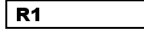
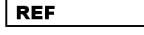
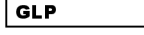
4. Tietz NW, ed. Clinical Guide to Laboratory Tests. 3rd ed. Philadelphia: WB Saunders 1995:919.
5. Tietz Fundamentals of Clinical Chemistry. 5th ed. Burtis CA, Ashwood ER, eds. Philadelphia, PA: WB Saunders Company; 2001:605.
6. Clinical and Laboratory Standards Institute (CLSI). Evaluation of Stability of In Vitro Diagnostic Reagents; Approved Guideline. CLSI Document EP25-A. Wayne, PA: CLSI; 2009.
7. Clinical and Laboratory Standards Institute (CLSI). Defining, Establishing and Verifying Reference Intervals in the Clinical Laboratory; Approved Guideline – Third Edition. CLSI Document EP28-A3c. Wayne, PA: CLSI; 2010.
8. Clinical and Laboratory Standards Institute (CLSI). Protocols for Determination of Limits of Detection and Limits of Quantitation; Approved Guideline. CLSI Document EP17-A. Wayne, PA: CLSI; Vol. 24 No. 34.
9. Clinical and Laboratory Standards Institute (CLSI). Evaluation of Precision of Quantitative Measurement Procedures; Approved Guideline – Third Edition. CLSI Document EP05-A3. Wayne, PA: CLSI; 2014
10. Passing-Bablok W et al. A General Regression Procedure for Method Transformation. J Clin Chem Clin Biochem 1988;26:783-79.
11. Clinical and Laboratory Standards Institute (CLSI). Method Comparison and Bias Estimation Using Patient Samples; Approved Guideline—Second Edition; Approved Guideline. CLSI Document EP09-A2. Wayne, PA: CLSI; Vol. 22 No. 19.
12. Clinical and Laboratory Standards Institute (CLSI). Interference Testing in Clinical Chemistry; Approved Guideline. CLSI Document EP07. Wayne, PA: CLSI; 3rd Edition. CHERIAN G., SOLDIN ST. Clin. Chem. 27/5:748-752 (1981)
13. Tietz Textbook of Clinical Chemistry, Second Edition, Burtis-Ashwood (1994).
14. Young DS. Effects of Drugs on Clinical Laboratory Tests. 3rd ed. Washington: AACC Press (1990).
15. Levinson S.S., Direct determination of serum chloride with a semiautomated discrete analyzer, Clin.Chem. 22:273-274, 1976
16. Wu, Alan H.B. Tietz Clinical Guide to Laboratory Tests. 4th ed. Saunders Elsevier, St. Louis, MO: 2006, 234-241.
17. Eisenman G. Glass Electrodes for Hydrogen and Other Cations, Principles and Practice. New York: Marcel Dekker Inc.; 1967:2.
18. Clinical and Laboratory Standards Institute (formerly NCCLS). Evaluation of Precision Performance of Quantitative Measurement Methods; Approved Guideline - Second Edition. Wayne, PA: Clinical and Laboratory Standards Institute; 2004. NCCLS Document EP05-A2.

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