

# MICROALBUMIN CALIBRATOR

# En

REF: VT-017    Microalbumin Calibrator    5 x 1 mL

*Changes made in the instructions for use are marked as grey.*

Package insert instructions must be carefully followed. Reliability of assay results cannot be guaranteed if there are any deviations from the instructions in this package.

### INTENDED USE

Validity Microalbumin Calibrator is for calibration of the Microalbumin Liquid assay.

### CONTENTS / MATERIALS PROVIDED

Microalbumin Calibrator  
**REF: VT-017**  
 Package: 5 x 1 mL **Lyophilized**

### For use with:

Validity Microalbumin Liquid Reagent  
 Validity Microalbumin Control Level I  
 Lyophilized  
 Validity Microalbumin Control Level II  
 Lyophilized

VT-017 contain the human serum.

Sodium Azide (0,09 %) is added as preservative.

### Materials Required But Not Provided:

1. Class A volumetric pipette for liquid transfer
2. Distilled or deionized water meeting the specifications equivalent to USP (United States Pharmacopeial Convention) purified water.

### STANDARDIZATION

The calibrator values were obtained by using the Validity Microalbumin reagent and Certified Reference Material (CRM) 470 (CRM 470) as the primary standard.

### CALIBRATOR STABILITY

Temperature-Conditions	Stability
Unopened at +2/+8°C	Expiry date on the vial.
Frozen and stored at -20°C	30 days

### LIMITATIONS

1. Please open the vial caps carefully. When you open, be careful not to scatter any powdery substance around or to escape from the vial.
2. Dissolve with distilled water with volume stated on the vial. Injector should not be used for the transfer process since there may be errors between 5-20% in liquid transfer with the injector. Use calibrated micropipettes.
3. Temperature of dry serum in the vial and distilled water must be +20/+25°C. After adding distilled water, close the vial cap tightly and store at +25°C around 5-10 minutes.
4. Wait for 30 minutes for dissolving process and mix thoroughly by gently inverting the vial at regular intervals, do not shake. Avoid formation of bubbles or foam. Protect from light. It is recommended to use a rotational mixer for routine mixing procedures.
5. After reconstitution, the calibrator serums are usually divided into small quantities (150-250 microliters) into Eppendorf tubes or sample cups of the device and stored in the refrigerator for freezing process. For serums prepared in this way, it is absolutely necessary to leave the serum at +25°C for 30 minutes before dividing it into small quantities. Do not refreeze after the serum is frozen and thawed once.
6. Calibrator serum precipitation is faster than normal serum. In order for the first and last parts to be homogeneous and to avoid precipitation, perform the process as fast as possible during separation.
7. The quality of the distilled water to be used in the dilution of the calibrator serum is very important. There may be significant deviations in the values due to bacterial contamination.
8. It is necessary to be careful against infectious agents in calibrator serum measurements.

### PREPARATION OF CALIBRATOR

Lyophilized serum calibrator should be reconstituted by adding distilled or deionized water with the amount stated on the label. Close the vial and wait for 30 minutes. Dissolve the contents of the vial by swirling gently to avoid the formation of foam. Do not shake.

### INDICATIONS OF INSTABILITY OR DETERIORATION

Presence of extreme turbidity or microbial growth may indicate deterioration.

## MICROALBUMIN CALIBRATOR

# En

### PRECAUTIONS



Human source material. Treat as potentially infectious material. Each plasma donor used in the preparation of this product has been tested by an FDA-approved method and found negative for the presence of HIV 1/2 HBsAg, HCV, HIV-Ag antibodies. However, none of the known testing methods can offer complete assurance that the hepatitis B virus, Human Immunodeficiency Virus (HIV) or infectious agents are not present. All human-based products should be handled in accordance with Good Laboratory Practice (GLP) principles using appropriate precautions. Safety data sheets are available at [www.validity.com.tr](http://www.validity.com.tr) or you may contact your local representative.

### WARNINGS

IVD: For in Vitro Diagnostic use only.

Do not use expired reagents.

Reagents with two different lot numbers should not be interchanged.

For professional use.

Follow Good Laboratory Practice (GLP) guidelines.

Contains sodium azide.

CAUTION: Human source samples are processed with this product. All human source samples must be treated as potentially infectious materials and must be handled in accordance with OSHA standards.

### Danger

H317 :May cause allergic skin reaction.

### Precaution

P280 :Use protective gloves / clothes / glasses / mask.

P264 :Wash your hands properly after using.

P272 :Contaminated work clothes should not be allowed to be used outside of the workplace.

### Intervention

P302+P352 :Wash with plenty of water and soap if it contacts with skin.

P333+P313 :Seek medical help if it irritates your skin or develops rash.

P362+P364 :Remove contaminated clothes and wash properly before using.

### Disposal

P501 :Dispose the vials and contents according to the local regulations.

### REFERENCES

1. Council Directive (2000/54/EC). Official Journal of the European Communities No. L262 from Oct. 17, 2000.
2. EU-Dir 1999/11 Commission Directive of 8 March 1999 adapting to technical progress the principles of Good Laboratory Practice as specified in Council Directive 87/18/EEC
3. Clinical and Laboratory Standards Institute, H26-A2, Validation, verification, and quality assurance of automated hematology analyzers; Approved Guideline - Second Edition.
4. Gabbay, K.H., Hasty, K., Breslow, J.L., Ellison, R.C., Bunn, H.F., and Gallop, P.M., J. Clin. Endocrinol. Metab. 44, 859 (1977).
5. US Department of Labor, Occupational Safety and Health Administration, 29 CFR Part 1910.1030, Occupational Exposure to Bloodborne Pathogens.
6. Ledue TB, Johnson AM. Commutability of serum protein values: persisting bias among manufacturers using values assigned from the certified reference material 470 (CRM 470) in the United States. Clin Chem Lab Med 2001;39(11):1129–1133.





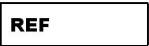










**Archem Sağlık Sanayi ve Tic. A.Ş. (With official contract based manufacturing agreement with Validity Sağlık Hiz. Sanayi)**  
Mahmutbey Mah. Halkalı Cad. No:124 Kat:4  
Bağcılar/Istanbul/Turkey  
Tel: + 90 212 444 08 92  
Fax: +90 212 629 98 89  
info@archem.com.tr www.archem.com.tr  
info@validity.com.tr www.validity.com.tr



# MICROALBUMIN CALIBRATOR

# En

SYMBOLS	
	In Vitro Diagnostic Medical Device
	Lot Number
	Calibrator
	Global Trade Item Number
	Reference Number
	Identifies Products to Be Used Together
	Product of Turkey
	Manufacturer
	Expiration Date
	Temperature Limits
	Consult Instructions for Use
	Caution
	Number of Tests

Validity